

# Your First IQCP Inspection



Edwina Szelag, MHM, BSMT(ASCP)

May 13, 2016



# Your Speaker

**Edwina Szelag MHM, BSMT(ASCP)**  
Technical Consultant II, Florida

- Board Certified ASCP
- Co-Chair and Co Founder Central Florida POC Network
- POCC of the Year 2006
- Clinical Laboratory Consultant





# Agenda

**1** Updates on IQCP from CAP, CMS, COLA, and TJC

**2** Observations from Completed IQCP Inspections

**3** IQCP Inspection Helpful Hints

**4** IQCP Inspection Question and Answer Session



# Updates on IQCP from CAP, CMS, COLA, and TJC



# CAP IQCP Resources

## IQCP Resources FAQ (56 questions)

- Check webpage for updates.

## POC Checklist 2015

## All Common Checklist 2015

- Read every line and have all the documentation listed under notes and comments.

## Instructions for Inspecting IQCP

## 2015 IQCP Requirements

## IQCP Inspector Tip Sheet

## IQCP Dos and Don'ts



# Highlights of CAP IQCP Resources FAQs

If my instrument has a control process that uses liquid control materials on-board the instrument or within a test cartridge, or uses a device, such as an optical filter or electronic control simulator, do I need to implement an IQCP to meet daily QC requirements?  
(UPDATED 11/23/2015)

- The default CLIA regulations were **written for the traditional daily testing of two levels of external control materials**. To be considered an external control material, the control material must follow the entire testing process, from **sample introduction through the analytic pathway**. It must also be a different type of material or from a different lot number than used to calibrate the instrument. Laboratories must carefully evaluate the control processes used to determine if they control the full analytic testing process. If the control process does not meet the criteria described for external control materials, the laboratory must either perform additional QC testing using appropriate external control materials or implement an IQCP to meet daily QC requirements.

Data sourced from CAP IQCP FAQs Dec. 18, 2015. Permission granted by  
College of American Pathologists.



# Highlights of CAP IQCP Resources FAQs

What are the QC requirements if I have multiple identical instruments/devices/cartridges in use but do not wish to develop an IQCP? (UPDATED 11/23/2015)

- Generally, the CAP and CLIA require at least two levels of external **QC for each device and cartridge, each day of testing**. Different CAP and CLIA requirements exist in some discipline and subsdiscipline areas (eg, coagulation, blood gases, and microbiology). The QC requirements for nonwaived testing, as written in the 2015 checklist edition, must be followed if an IQCP is not implemented by January 1, 2016.

Is a separate Risk Assessment required for each site if the same instrument/device/test is used in multiple areas within a CAP number?

- No. The laboratory has an option. Individual assessments may be performed or a single risk assessment (RA) may be used when there are multiple sites performing testing under a single CAP number. If a single RA is performed, all variations in the required components must be taken into account when conducting the RA (eg, **differences in sites, environments, or personnel**). A laboratory can then develop one IQCP that accounts for all of the differences in the RA or can develop individual IQCPs to address differences by site. Each device used must be monitored in some way, as well as each location.

Data sourced from CAP IQCP FAQs Dec.15, 2015. Permission granted by College of American Pathologists.



# Highlights of CAP IQCP Resources FAQs

My instrument (or kit) manufacturer provides risk assessment information for implementing an IQCP. Is this acceptable to use?

- Yes. It is acceptable to use information provided by an instrument or kit manufacturer as a supplement in the risk assessment.
- The laboratory needs to perform its own evaluation of all five elements of risk.
- **Laboratories cannot use manufacturer risk assessment alone**

Data sourced from CAP IQCP FAQs Dec. 18, 2015. Permission granted by College of American Pathologists.





# Highlights of CAP IQCP Resources FAQs

For the external QC at least every 31 days, how many levels of controls need to be run? If multiple devices are in use, can we run the external QC using a subset of devices?

- The laboratory must define the control procedures to be followed based on the risk assessment performed.
- The decision on the number of controls needed and the **use of subsets of devices** using the same reagent lot if multiple devices are used may be defined by the laboratory in the quality control plan, if appropriate, and be approved by the laboratory director based on the risk assessment evaluation and the supporting data used in the risk assessment

Data sourced from CAP IQCP FAQs Dec. 18, 2015. Permission granted by College of American Pathologists.



# Highlights of CAP IQCP Resources FAQs

Will the checklist requirements for more frequent QC for some types of testing still apply (eg, coagulation, blood gases) if a laboratory implements an IQCP?

- During the Risk Assessment process for a test that is eligible for IQCP, the laboratory must evaluate the potential sources of errors, manufacturer's instructions, and historical test performance to identify the appropriate control processes. The laboratory's Quality Control Plan **may define a frequency less than the minimum frequency defined in the CAP checklist** if it is determined to be acceptable based on the risk assessment.
- **If an IQCP is not implemented, the minimum QC frequency defined in the CAP checklists and default CLIA requirements must be followed.**
- In all cases, manufacturer's requirements for QC must be followed, at a minimum.



# Highlights of CAP IQCP Resources FAQs

What are the QC requirements if I have multiple identical instruments/devices/cartridges in use but do not wish to develop an IQCP?

- Without an IQCP, the existing CLIA and 2015 CAP checklist requirements will apply; generally, at least **two levels of external QC for each device and cartridge**, each day of testing (or more frequently as specified in a discipline or subspecialty).

What is required for ongoing assessment of an IQCP?

- Ongoing assessment must include evaluation of errors relating to the different phases of the testing process, QC failures and corrective action, and complaints from clinicians and other providers on the quality of results. It must also include a determination of the need to reassess and revise the IQCP.
- Quality control and instrument/equipment maintenance and function check data must continue to be reviewed at least monthly.
- Additionally, each IQCP must be **assessed annually** for effectiveness and revised, as necessary.



## Highlights of CAP IQCP Resources FAQs

Does the quality assessment monitoring for IQCP need to be included in the quality management (QM) program?

- If used, IQCP must be incorporated into the quality management program. Ongoing quality assessment of an IQCP must include evaluation of errors relating to the different phases of the testing process, QC failures and corrective action, complaints from clinicians and other providers on the quality of results, and an **annual assessment** of the effectiveness of the IQCP. Some of these items are often included in the QM plan already. The laboratory may consider including ongoing assessment of **these items as quality indicators.**

Data sourced from CAP IQCP FAQs Dec. 18, 2015 Permission granted by  
College of American Pathologists.



# Highlights of CAP IQCP Resources FAQs

## How will IQCP be inspected?

- Inspectors will look for compliance with the requirements defined in the 2015 checklist for IQCP.
- For each CAP number, requirements for compliance will include:
- Risk Assessment, including evaluation of all of the following:
  - All five required elements (Reagents, Environment, Specimen, Test System, Testing Personnel)
  - All phases of testing: pre-analytic, analytic, and post-analytic
  - Data from the laboratory's own environment, instrument/equipment performance, and testing personnel
  - All variations in test performance (eg, multiple test sites, devices, types of testing personnel, etc)



# Highlights of CAP IQCP Resources FAQs

## How will IQCP be inspected (continued)?

- Written Quality Control Plan defining types of control processes used, criteria for acceptable performance, and frequency evaluated. QC may not be performed less frequently than defined in the manufacturer's instructions.
- Approval of the written IQCP by the laboratory director prior to implementation (signed and dated)
- Ongoing assessment of errors, QC failures, and complaints, including the need to reassess the risk assessment and quality control plan
- **Annual review** of each IQCP
- **Use of CAP forms to maintain a list of Individualized Quality Control Plans and a summary of each IQCP**
- Inspectors **may cite deficiencies** when any of the above elements are not in compliance with checklist requirements. **The decision on whether the level of risk for any of the elements evaluated in the risk assessment is acceptable is left to the discretion of the laboratory director.**



# Read Every Line in the Requirement Text and Notes

## Example: CAP Checklist COM.50300 Risk Assessment (RA)

The QC study performed to assess the performance and stability of the tests must support the QC frequency and elements defined in the laboratory's quality control plan. The study must include data representing, **at a minimum the maximum interval between runs of external quality control**. The laboratory may use historical data during the risk assessment for tests already in place.



# How To Get Help!

From: Laboratory Accreditation Program-General [<mailto:accred@cap.org>]

Sent: Friday, February 19, 2016 3:27 PM

To: Smith, Jane; Smith, Jane

Subject: SR 2460382 Re: Interpretation of the Requirement Text COM.50300 [REF:43770134012]

Hello Jane,

The Lab Director should be involved in the design function of the IQCP study. They must review and sign the study before implementation. CAP is not restrictive of the tool used to perform the IQCP, however, it must include all five required elements (reagents, environment, specimen, test system, and testing personnel). All phases of testing must be addressed (pre-analytic, analytic, and post-analytic). Data from the lab's own environment, instrument/equipment performance, and testing personnel must be used. All variations in test performance, for example multiple test sites, devices, types of testing and personnel must be studied. Concerning your question about whether CAP has a recommended study for QC, there are so many variables involved with the labs that are CAP accredited, each site must develop their own protocol to assess the risk at their facility.

Thank you for your question.


Sincerely,

Senior Inspection Specialist






# IQCP Inspector Tip Sheet

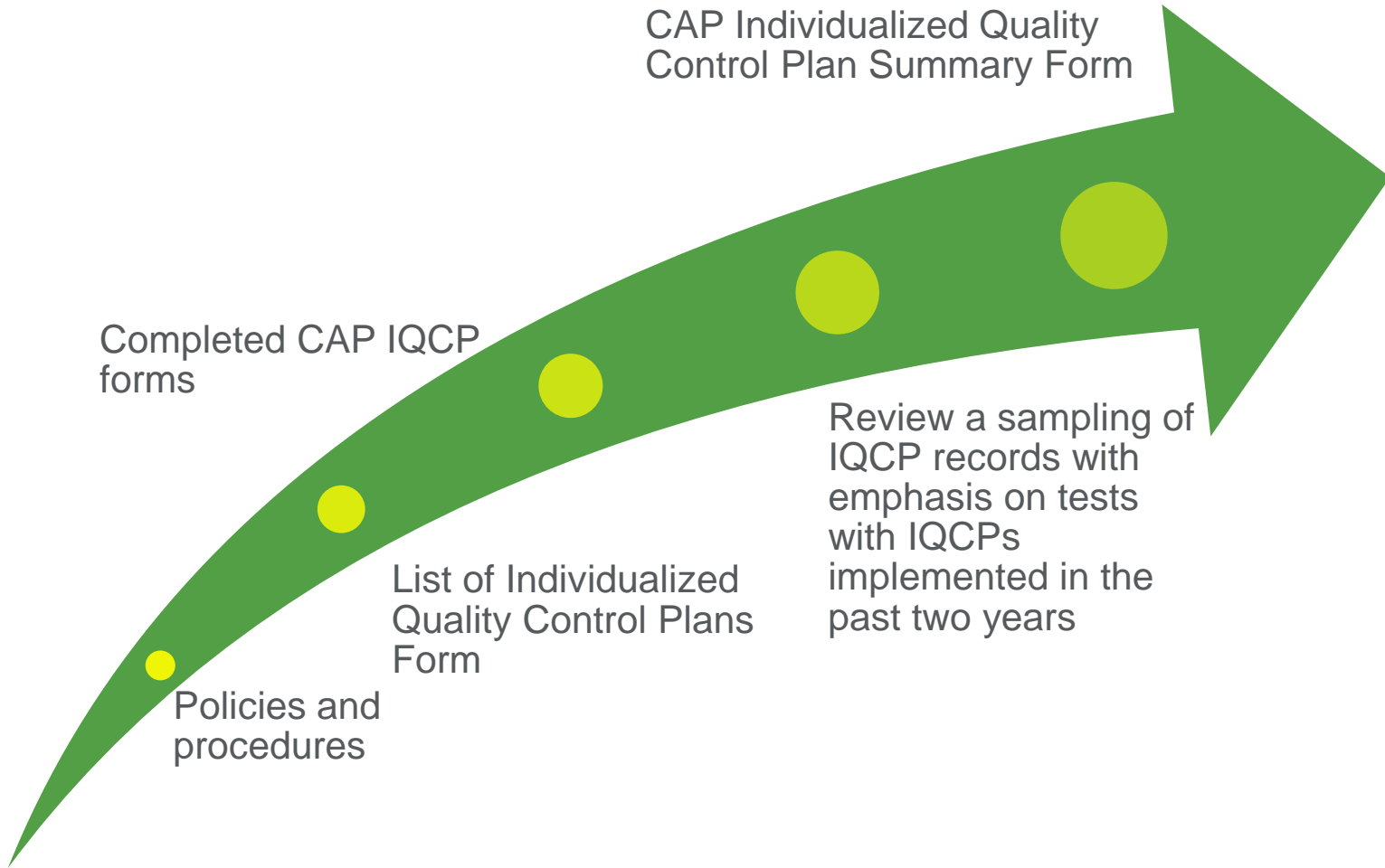

Data sourced from IQCP Inspector Tip Sheet Sept. 2015. Permission granted by College of American Pathologists.



# IQCP Inspector Tip Sheet for Risk Assessment (RA)

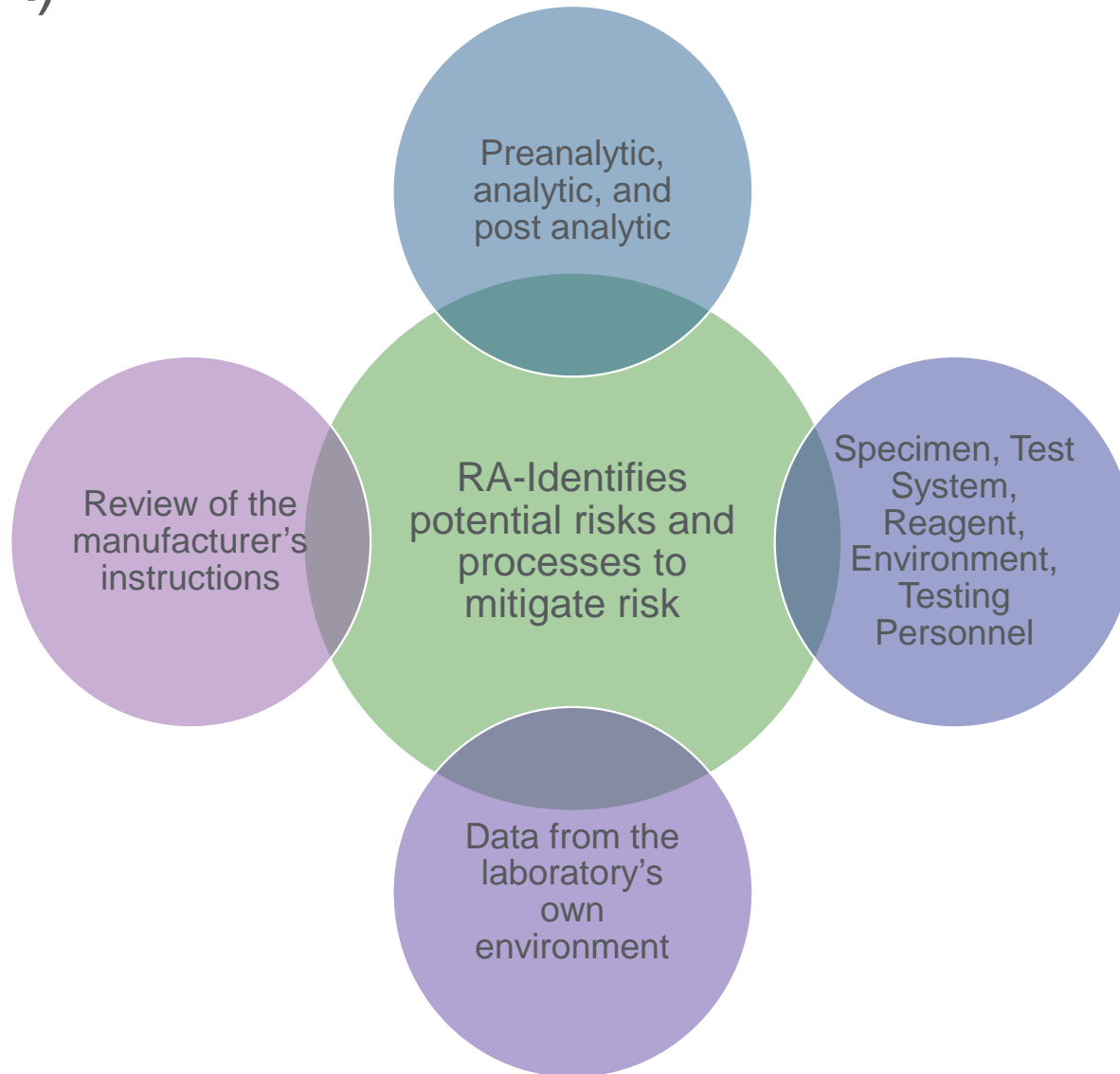


Data sourced from IQCP Inspector Tip Sheet Sept. 2015. Permission granted by College of American Pathologists.

FOR EXTERNAL USE PRINT AND DISTRIBUTION ARE PERMITTED



# Inspector Tip Sheet for Quality Control Plan (QCP)

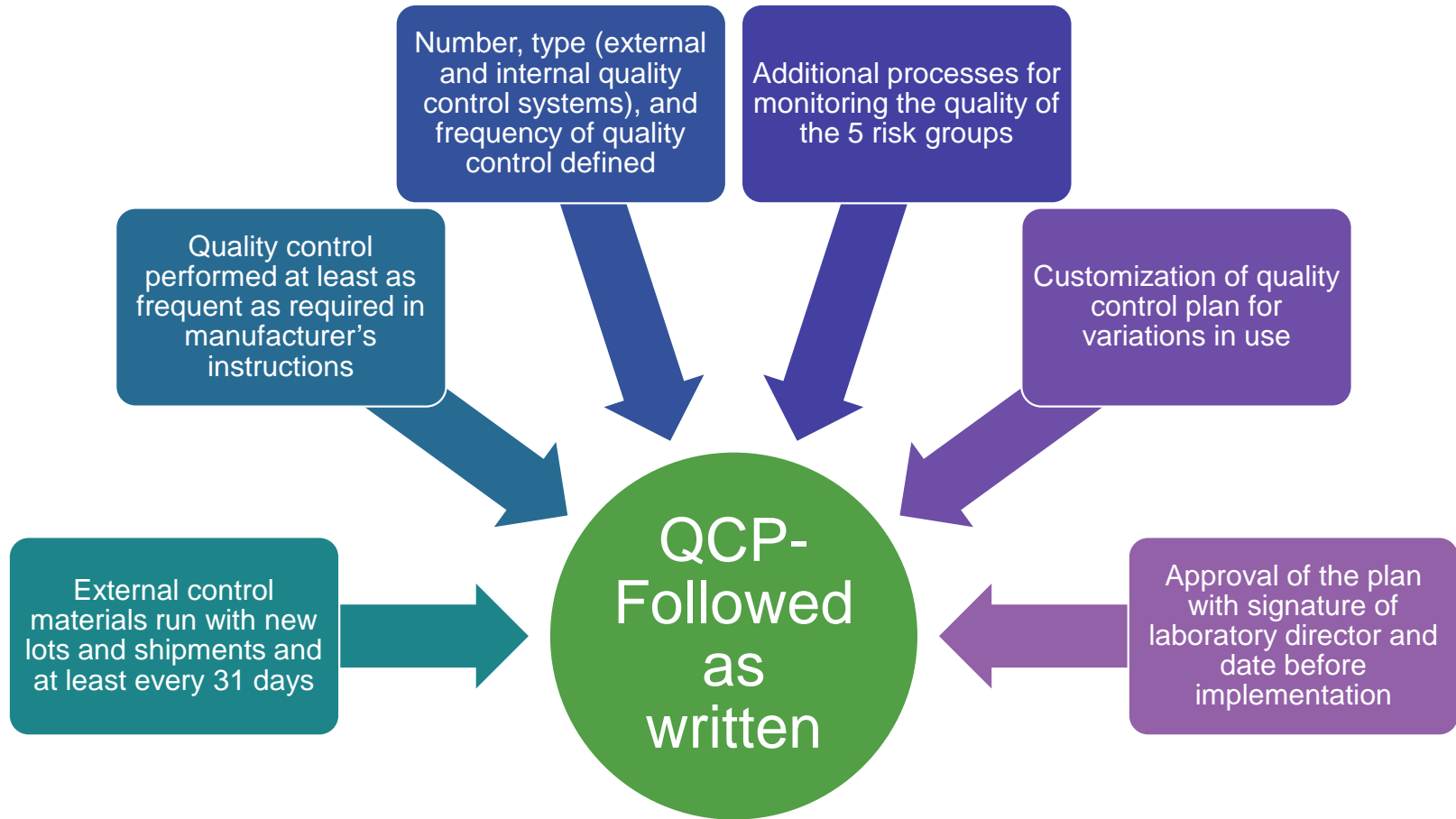

Data sourced from IQCP Inspector Tip Sheet 2015. Permission granted by College of American Pathologists



# IQCP Inspector Tip Sheet for Quality Assessment (QA)



Data sourced from IQCP Inspector Tip Sheet Sept 2015. Permission granted by College of American Pathologists.



## Inspector IQCP Do's and Don'ts

IQCP REQUIREMENT	DO CITE IF:	DON'T CITE BECAUSE:
COM.50300	1.) Risk Assessment (RA) is missing one or more of the five required components (specimen, reagent, environment, testing personnel, test system)	1.) The format of RA is not "user-friendly" - RECOMMEND
	2.) RA doesn't cover all three phases of testing: pre-analytic, analytic, and post-analytic	2.) The RA doesn't look like the ones in YOUR lab - DISCUSS
	3.) RA did not include in-house data (previous QC records, environmental monitoring, etc.) or did not involve laboratory personnel	3.) You disagree with the acceptability of a specific risk - DISCUSS

Permission granted to use Inspector IQCP Do's and Don'ts Sept. 2015 Attachment E by College of American Pathologists.



## Inspector IQCP Do's and Don'ts

<b>IQCP REQUIREMENT</b>	<b>DO CITE IF:</b>	<b>DON'T CITE BECAUSE:</b>
COM.50400, COM.50500	4.) Quality Control Plan (QCP) was not signed by the laboratory director prior to implementation	4.) You disagree with the frequency of the QC being run - RECOMMEND
	5.) QC is performed less frequently than specified in manufacturer's instructions	5.) You think the QCP does not address potential risks - RECOMMEND
	6.) QCP is not followed as written	6.) You disagree with the acceptability of QC to mitigate a specific risk - DISCUSS

Permission granted to use Inspector IQCP Do's and Don'ts Sept. 2015 Attachment E by College of American Pathologists.



## Inspector IQCP Do's and Don'ts

IQCP REQUIREMENT	DO CITE IF:	DON'T CITE BECAUSE:
COM.50600	7.) Quality Assurance process does not monitor devices used in all locations	7.) You don't think that the lab has adequately addressed potential patient outcomes - RECOMMEND
	8.) Serious quality concerns or adverse patient outcomes have not been addressed – MAY ALSO NEED TO CITE TLC.10460	

Permission granted to use Inspector IQCP Do's and Don'ts Sept. 2015 Attachment E by College of American Pathologists.



COLLEGE of AMERICAN  
PATHOLOGISTS

# Individualized Quality Control Plan Summary

Complete a separate form for each IQCP in use and present to the inspector during the on-site inspection.

Laboratory Name:  Laboratory Section/Department:  CAP Number:

1) Instrument/Device Include name, manufacturer, and model	2) Tests List all tests included under the IQCP	3) Number of Devices In Use	4) List of Test Sites* If used in more than one area	Date of Director Approval	Date Implemented	Date Retired
				Click here to enter a date.	Click here to enter a date.	Click here to enter a date.

## 5) Control Processes Used to Monitor Risk

*Include a brief statement about each control process – list the monitor and frequency evaluated*

Reagents	Environment	Specimen	Test System	Testing Personnel	Other

This is not to be used for risk assessment data. Only quality control measures.


The form is intended to be used for developing an IQCP or the performance for a risk assessment. The form is to be completed by laboratories preparing for a CAP inspection. Inspectors will use this form and IQCP List as tools for audit the IQCPs in use during a CAP onsite inspection. Permission granted to use this form by College of American Pathologists.





Many Great Pages Worth Reading.

This is the most recent document available on CMS.gov IQCP webpage

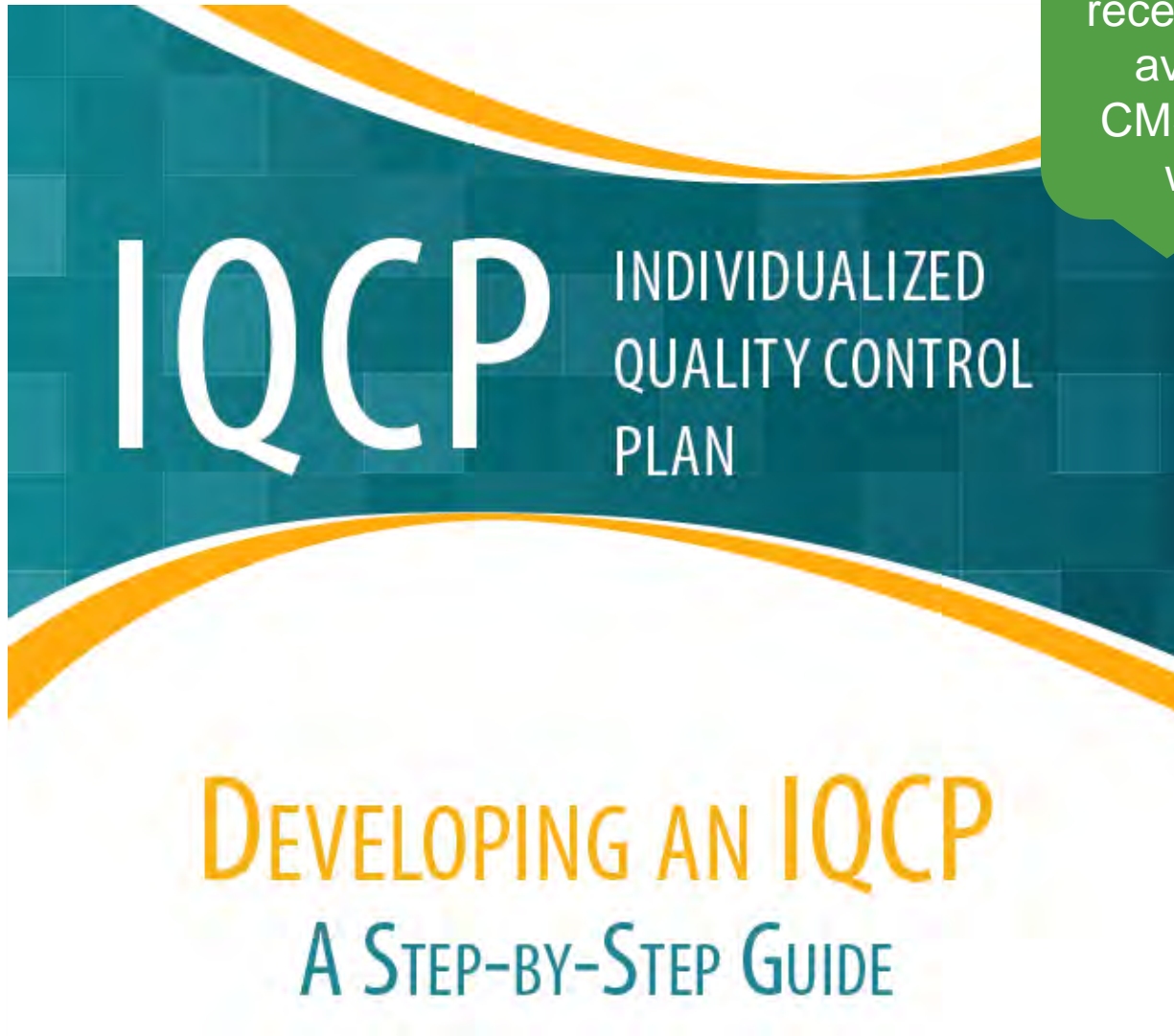

<http://www.cms.gov/Regulations-and-Guidance/Legislation/CLIA/Downloads/IQCP-Workbook.pdf>



# Quality Control Plan Tips to Remember


## A complete QCP must:

- Provide for immediate detection of errors for each phase of the testing process (i.e. before, during, and after testing) for the test.
- Specify the number, type, and frequency of testing QC material(s).
- Contain criteria to determine acceptable QC results.
- Require the laboratory perform QC as specified by the manufacturer's instructions, but not less than the manufacturer's instructions.
- Indicate that your Laboratory Director reviewed, signed, and dated the QCP document.
- **If your QCP does not address all five items listed above, you do not have a QCP.**
- **Go back and investigate what is missing.**

Data sourced from: <http://www.cms.gov/Regulations-and-Guidance/Legislation/CLIA/Downloads/IQCP-Workbook.pdf> pg 42



# Quality Assessments - Documents to Consider

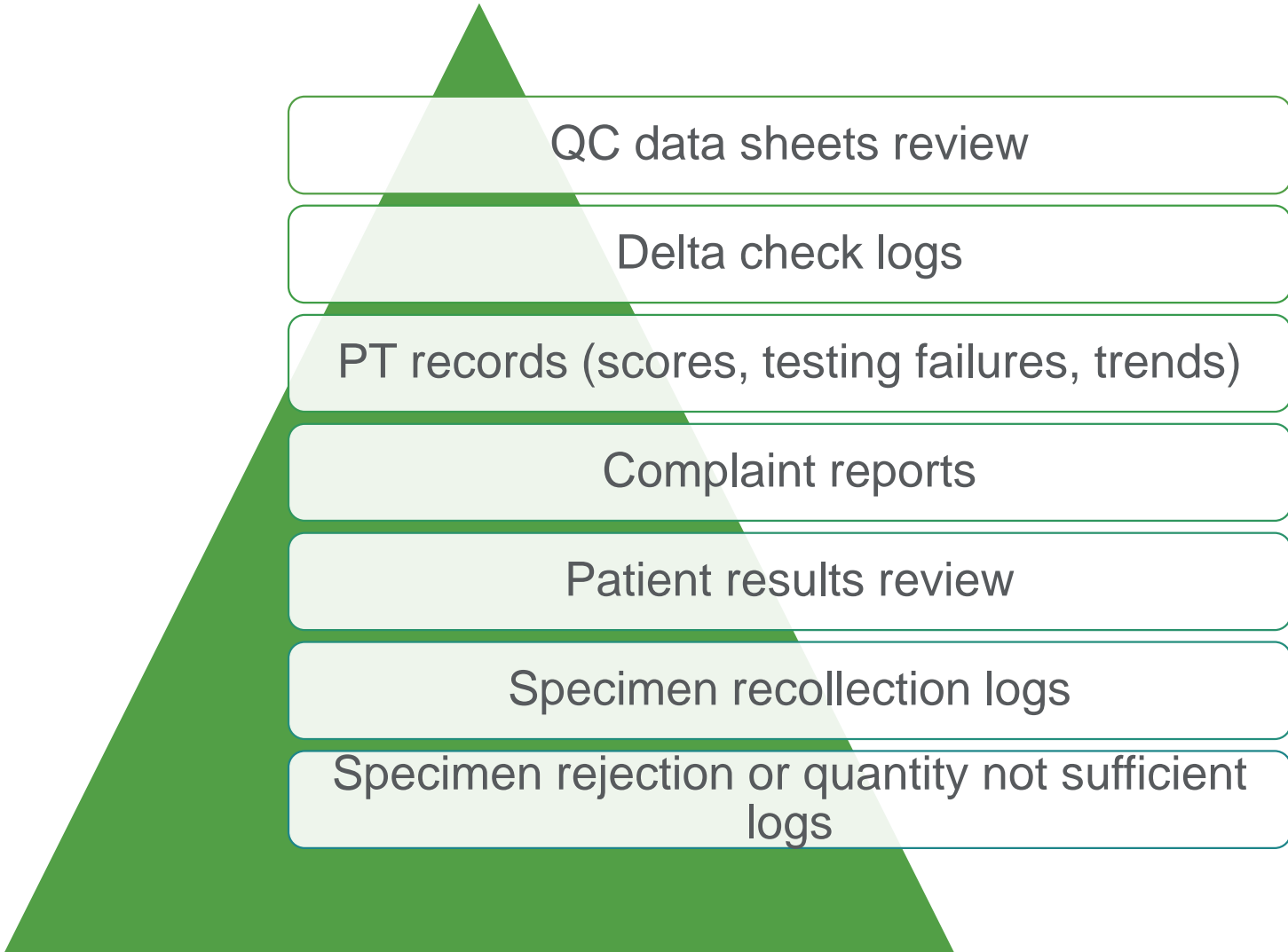

QC data sheets review

Delta check logs

PT records (scores, testing failures, trends)

Complaint reports

Patient results review


Specimen recollection logs

Specimen rejection or quantity not sufficient logs

Data sourced from <http://www.cms.gov/Regulations-and-Guidance/Legislation/CLIA/Downloads/IQCP-Workbook.pdf> pg 45





# Quality Assessments - Documents to Consider

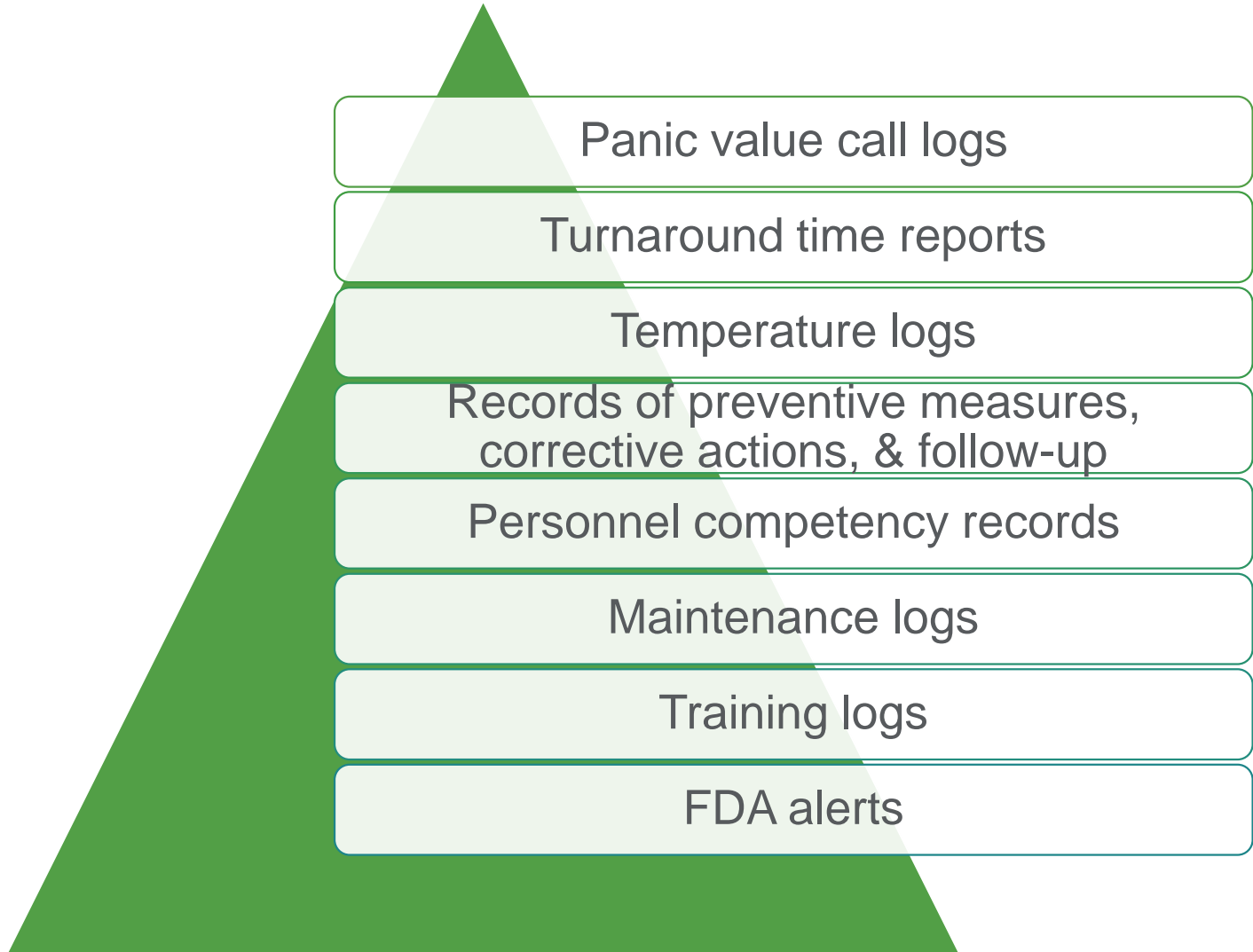

Data sourced from <http://www.cms.gov/Regulations-and-Guidance/Legislation/CLIA/Downloads/IQCP-Workbook.pdf> pg 45

# Quality Assessment (QA) Helps You.....

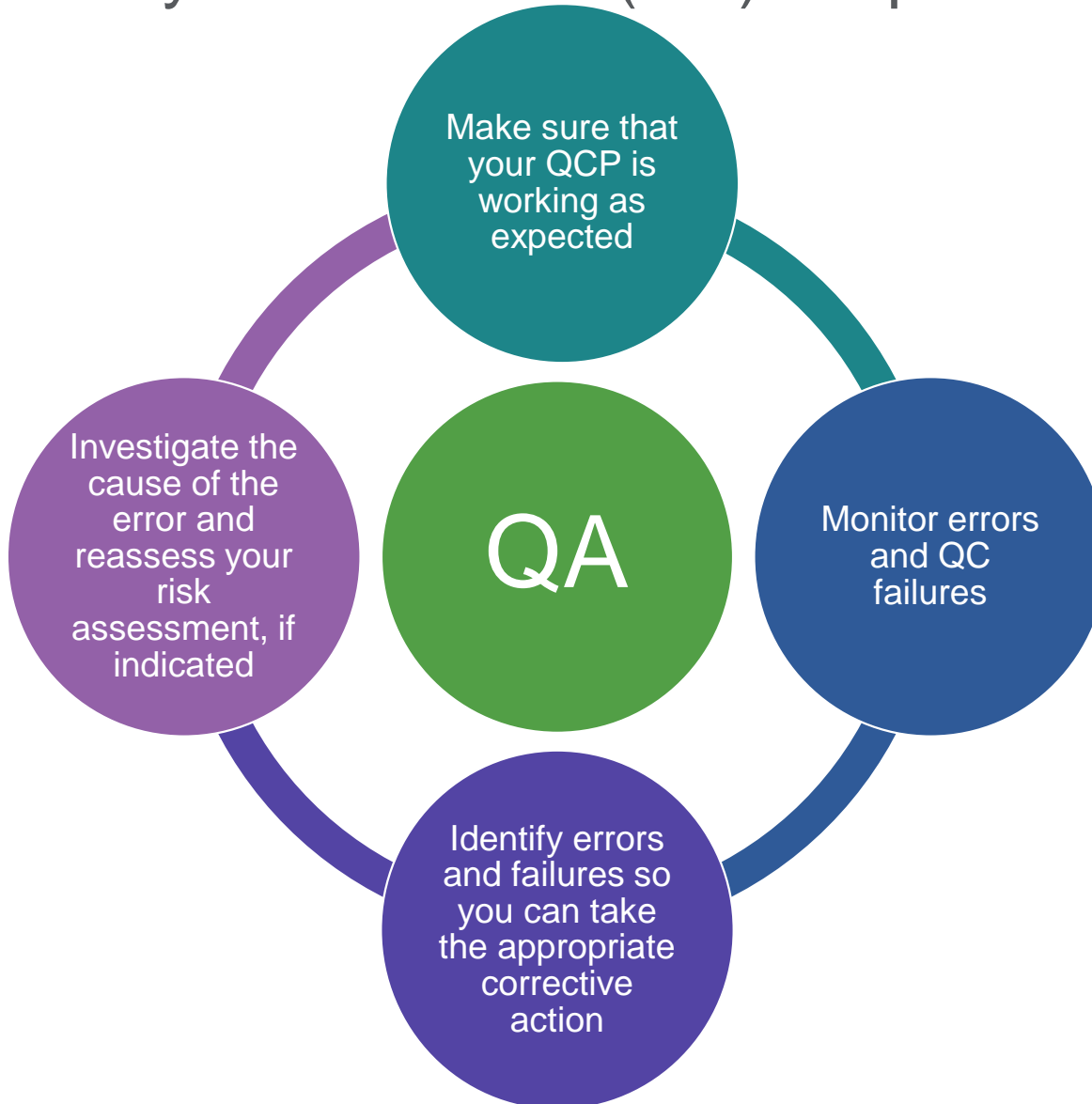

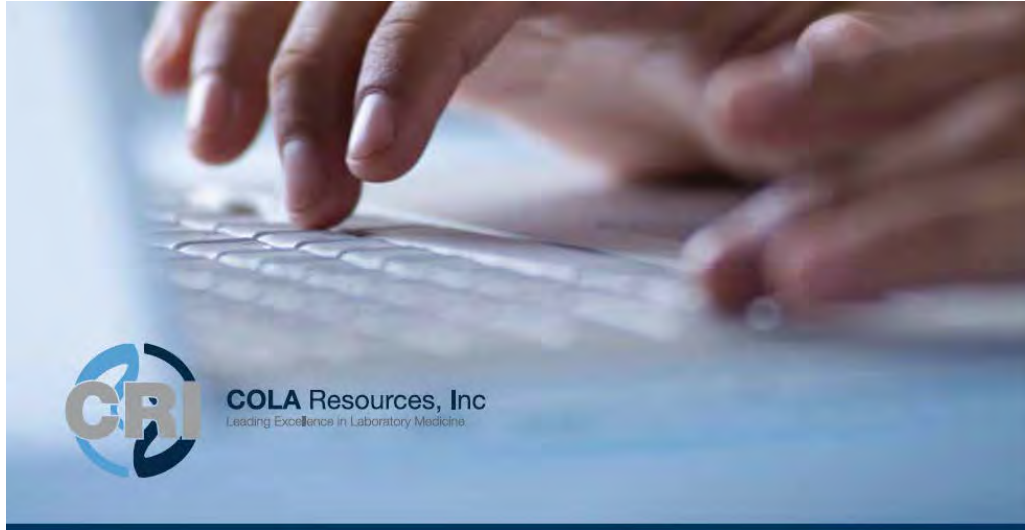

# Individualized Quality Control Plan (IQCP)

IMPLEMENTATION GUIDE



Permission to Use Copyrighted Material granted by:  
COLA Resources, Inc., (CRI) a laboratory medicine educational and  
consulting organization, located at 9881 Broken Land Parkway, Suite 215  
Columbia Maryland 21046



# The Joint Commission (TJC) Update IQCP is Not Required

	Jan. 1, 2016
Default (2-3 levels external QC/day)	<input checked="" type="checkbox"/>
EQC (Equivalent QC)	<input type="checkbox"/>
IQCP	<input checked="" type="checkbox"/>



# Observations from Completed IQCP Inspections



# CLIA Inspector Comments

## Inspection Observations:


- I've only seen a few IQCPs implemented, but they were good. I had a couple of proactive labs and they were prepared.

## Know What Tests Qualify for IQCP

### Inspection Observations:

- Labs are unfamiliar with the eligibility determination requirement.
- The lab had incorporated or listed waived tests on their form which are ineligible tests

### Eligibility Determination for Individualized Quality Control Plan (IQCP) Option

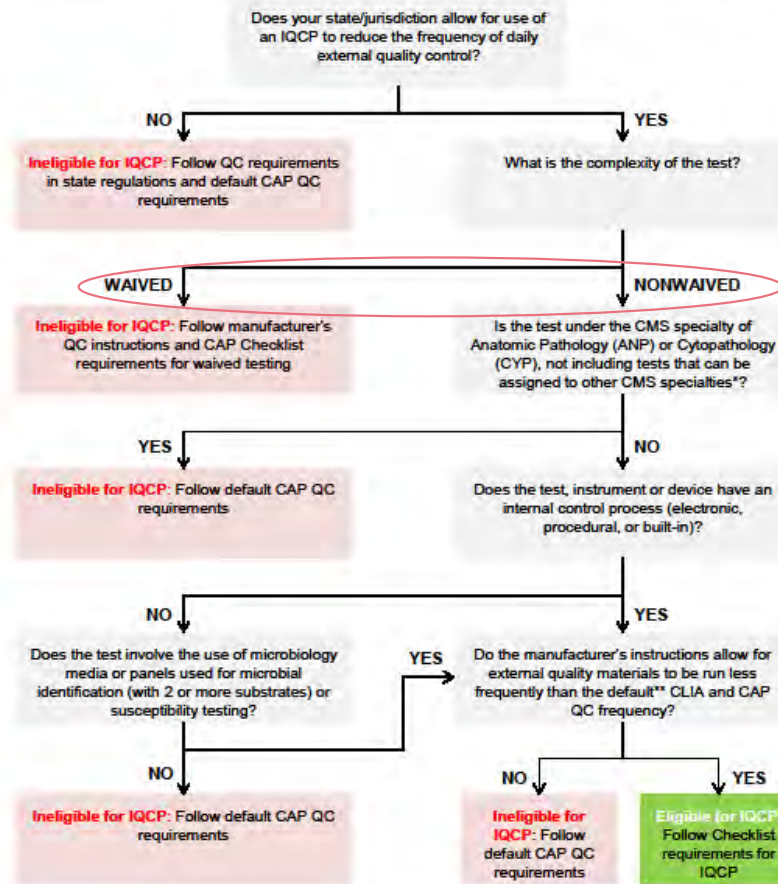


\* ANP or CYP tests are ineligible for IQCP unless the testing can be billed under another CMS specialty.


\*\* The default CAP QC frequency for external quality control materials is as follows:

- 1) Quantitative tests - two controls at different concentrations each day of patient testing, except for Coagulation tests (two levels every eight hours) and Blood Gas testing (one level every eight hours)
- 2) Qualitative tests - positive and negative controls each day of patient testing.



### Eligibility Determination for Individualized Quality Control Plan (IQCP) Option

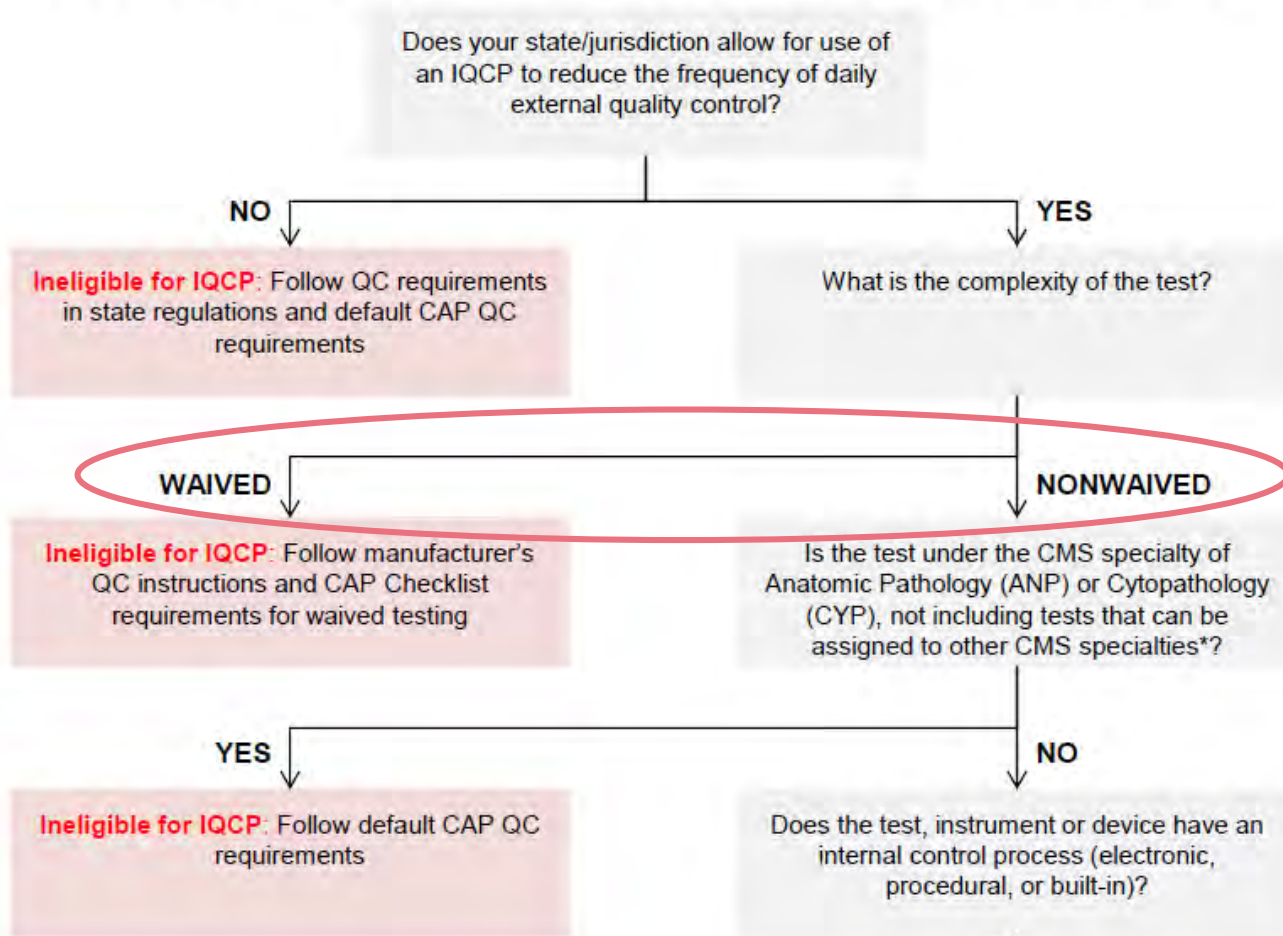

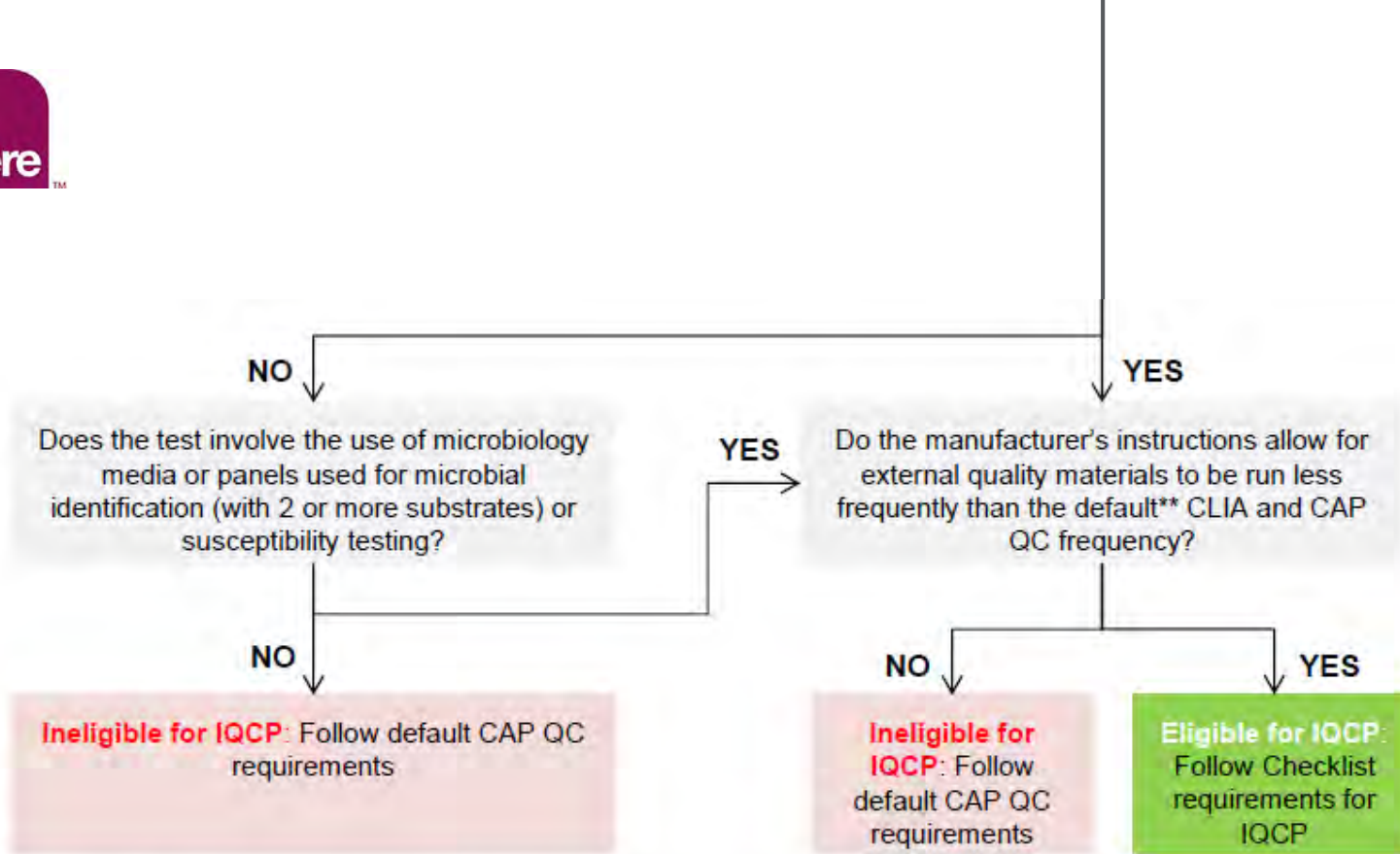

\* ANP or CYP tests are ineligible for IQCP unless the testing can be billed under another CMS specialty.

\*\* The default CAP QC frequency for external quality control materials is as follows:

- 1) Quantitative tests - two controls at different concentrations each day of patient testing, except for Coagulation tests (two levels every eight hours) and Blood Gas testing (one level every eight hours)
- 2) Qualitative tests – positive and negative controls each day of patient testing.



# CAP Checklist COM.50200 IQCP Test List/Summary

Individualized Quality Control Plan Summary

List of Individualized Quality Control Plans

- Inspection Observations:

These documents should be available as stand alone documents for the inspection team.

In the summary plan, test systems were not listed.



# CAP Checklist COM.50300 Risk Assessment (RA)

## Inspection Observations:

- Laboratories have received IQCP support documents from manufacturers. They have not added their own lab's data and included this information in their plan.
- Laboratories have not evaluated potential sources of error in preanalytic, analytic and postanalytic testing process.



# CAP Checklist COM.50300 Risk Assessment (RA)

## Inspection Observations:

- Laboratory was too vague on addressing potential risks associated with the instrument being used in multiple environments. The RA with multiple identical devices must show an evaluation was performed if there are differences in testing personnel and or environments where testing is performed.





# CAP Checklist COM.50300 Risk Assessment (RA)

## Inspection Observations:

- POCC was only person involved in writing IQCP. The laboratory must involve a representative sample of testing personnel in the process of conducting the RA.



# CAP Checklist COM.50300 Risk Assessment (RA)

## Inspection Observations:

- Laboratory's risk score was high for a risk, and when asked what was done for mitigation they didn't have an answer. Laboratory must have documentation for what do you do to reduce mitigating factor



# CAP Checklist COM.50500 Quality Control Elements

## Inspection Observations:

- Frequency and type of QC has not been transferred from the test procedure and specified in their IQCP.



# CAP Checklist COM.50500 Quality Control Elements

## Inspection Observations:

- IQCP does not include manufacturer's instructions to ensure the frequency of QC is not less than required by manufacturer. Package inserts should be available as part of the laboratory's documentation to ensure frequency of QC is appropriate.



# CAP Checklist COM.50600 Quality Assessment Monitoring

## Inspection Observations:

- Some labs have not specified how IQCP will be incorporated into their QA programs, and how it will be evaluated/corrective actions if needed. Review and documentation of QC and corrective actions is being done but documentation in QA monitoring has been incomplete.

# Laboratory Comments Post Inspection

## Laboratory Comments about IQCP Inspections

- The IQCP was good enough for our laboratory director to sign how was it lacking? Laboratory thought that was the purpose of the "individual" in the IQCP.
- The inspector was very knowledgeable and respectful. The inspector was very patient and even provided suggestions on how to improve my IQCP.
- We did not have to revert to doing 2 or 3 levels of QC daily or every 8 hours because we had a signed IQCP even though our inspector thought our IQCP was lacking in areas.



# Documentation Ideas.....

<b>Individual Quality Control Plans:</b>	<b>Document Retention Period</b>
<b>Risk assessments, validation data, approved quality plans and ongoing quality assessments. These documents must be retained for two years following the discontinuing of testing or the Individual Quality Control Plan.</b>	<b>2 years</b>



# IQCP Inspection Helpful Hints





# Top Tips

Read all the CAP  
FAQs

CAP FAQs are  
updated frequently

Use Email Contacts

Best Guide:  
Developing an IQCP  
A Step-by-Step  
Guide

Read CAP  
Inspecting the New  
IQCP

Training,  
Competency, and  
Educational  
Documentation

Read Every Line of  
the CAP Checklist  
Including  
Requirement Text  
and Notes

Extra External  
Quality Control  
Materials

Read AACC POC  
List Serve

Watch EALERTS  
from CAP

E-Alerts from TJC



## For More Information

CMS	<a href="mailto:IQCP@cms.hhs.gov"><u>IQCP@cms.hhs.gov</u></a>
CAP	<a href="mailto:accred@cap.org"><u>accred@cap.org</u></a>
Joint Commission	<a href="mailto:qualitylabs@jointcommission.org"><u>qualitylabs@jointcommission.org</u></a> Go to the Leading Practice Library
COLA	<a href="mailto:info@cola.org"><u>info@cola.org</u></a>
Alere	<a href="mailto:IQCP@alere.com"><u>IQCP@alere.com</u></a> <a href="http://www.alere.com/IQCP"><u>www.alere.com/IQCP</u></a>



# IQCP Inspection Question and Answer Session



9975 Summers Ridge Road  
San Diego, CA 92121  
**Tel:** 877.441.7440  
**Fax:** 877.441.7441  
**Web:** [www.alere.com](http://www.alere.com)

# Thank you!

© 2016 Alere. All rights reserved. The Alere Logo, Alere, and Knowing now matters are trademarks of the Alere group of companies. SAHIGHR-0070-B